
尺寸: W137mmH149mm

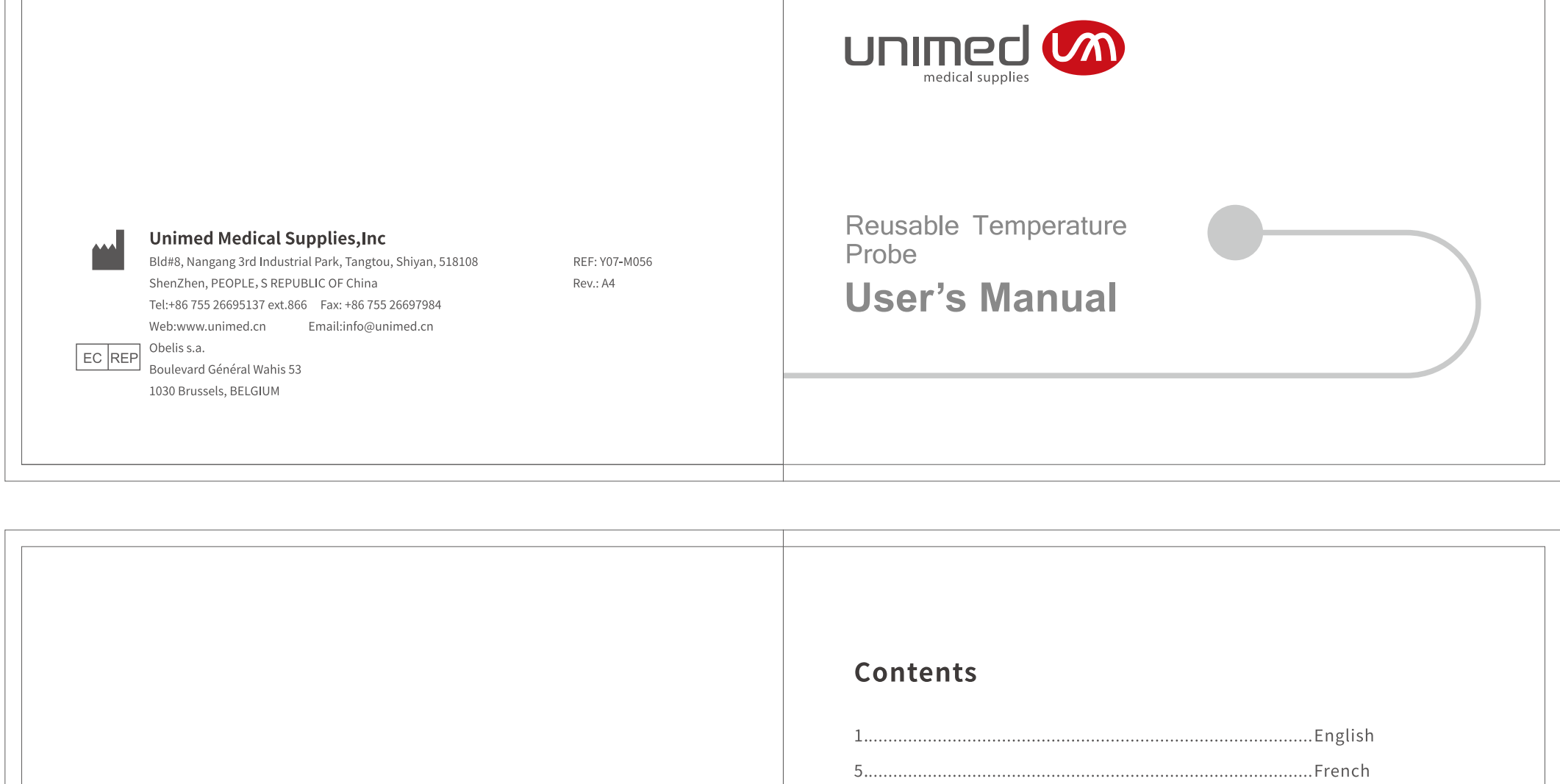unit or other RF source. If probes must be used together with electro-surgical apparatus, the instrument cables must be connected to the patient. The patient must be connected to the patient ground. The patient and the probes are connected should be checked for adequate isolation from electrical grounds at radio frequencies.

Caution:

Overheating of the probes could result in damage to internal wires and loss of electrical isolation or improper measurement readings.

Operation procedures:

To ensure the compatibility of the probe, see the attached the compatible label to ensure the device is compatible with the temperature probe.

Make sure the probe plug matches the connector on the device. Then fully plug the probe into the device. The probe is used to measure the temperature of the tissue. The probe is used to measure the temperature of the tissue. The probe is used to measure the temperature of the tissue.

Place the probe tip on the site where is intended to be measured. Fix the probe tip with medical tape or a bandage.

Read the measurement readings on the display.

Start to measure and get the measurement readings after the response time (see above).

Disinfections:
Probes should be cleaned before disinfections or sterilization to improve the effectiveness.


Probes may be disinfected by washing with 70% isopropyl ethanol. After washing, probes should be rinsed thoroughly with water. Brief immersion of the probe in detergent solutions is not harmful.

Never boil the T series reusable temperature probes.

Storage and handling:

When not in use, probes and leads should be loosely coiled and stored in room temperature. Don't wrap probes around equipment cases to avoid damaging internal wires.

Environment operating range:

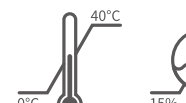Warranty:

Regular lifetime of Unimed temperature probes is two years. All T series reusable temperature probes carry one year warranty on workmanship, components and accuracy tolerances. Damage through misuse or modification is not covered. Probe life with normal use should exceed one year.

CAUTION:

Federal (U.S.) Law restricts this device to sale by or on the order of a physician.

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Description générale :

Symbol explanation

Cables	Lot number	Manufacturer	Manufacturer	Lotus free	Self-heating materials	Non-Chemical and Bio-Tissue Equivalents
						

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