
Intersorb™ Plus Soda Lime Spherasorb™ Soda Lime LoFloSorb™

Loose fill carbon dioxide absorbents for re-fillable absorbers

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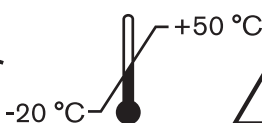en Instruction for use: Intersorb™ Plus, Spherasorb™ and LoFloSorb™

These instructions must be carefully read before using Intersorb Plus, Spherasorb and LoFloSorb. Also refer to the instructions of use provided by the anaesthetic machine manufacturer.

WARNING: A WARNING statement provides important information about a potentially hazardous situation which, if not avoided could result in death or serious injury.

CAUTION: A CAUTION statement provides important information about a potentially hazardous situation which, if not avoided could result in minor or moderate injury.

TECHNICAL DATA

Product Code	Absorbent Volume	Resistance to flow of 1 kg absorber at 60 lpm	Typical pH	Absorbent details
2172000	1.0 kg / 2.2 lb	0.7 mbar	Used below 10 Unused 12 – 14 depending on absorbent	For full details please contact Intersurgical
2173000	1.0 kg / 2.2 lb			
2174000	5 L			
2175000	5 L			
2178000	5 L			
2179000	5 L	0.8 mbar		
2180000	5 L			

INTENDED USE

To remove carbon dioxide from anaesthetic and respiratory gases delivered to a patient.

Intersorb Plus, Spherasorb and LoFloSorb Loose fill carbon dioxide absorbents are for the absorption of carbon dioxide within an anaesthetic breathing system and are compatible with all re-fillable absorbers designed for use with loose fill absorbents.

INTENDED USER

For use by qualified medical personnel.

CHOICES AVAILABLE

Product codes:

2172000 Spherasorb, pink to white colour change, 1kg bag

2173000 Spherasorb, white to violet colour change, 1kg bag

2174000 Spherasorb, pink to white colour change, 5L jerican

2175000 Spherasorb, white to violet colour change, 5L jerican

2178000 LoFloSorb, green to violet colour change, 5L jerican

2179000 Intersorb Plus, white to violet colour change, 5L jerican

2180000 Intersorb Plus, pink to white colour change, 5L jerican

For full details of Spherasorb, LoFloSorb and Intersorb Plus please contact Intersurgical or refer to the Intersurgical Absorbents catalogue.

PRINCIPLE OF OPERATION

These absorbents undergo a chemical reaction with the CO₂ within the breathing system. The hydroxides within the absorbent are converted to carbonates and water in an exothermic process. The CO₂ is therefore, bound within the absorbent.

WARNING If the moisture content of the absorbent is outside that set during the manufacturing process, this will affect the absorbing performance.

For more information on this issue please see the Safety Information section or contact Intersurgical.

FILLING RE-FILLABLE ABSORBERS

CAUTIONS

1. Fill the absorber in two or three stages. Fill to the filling mark of the absorber or if this is not present, to within about 5 mm of the top. Once the absorbent is filled into the absorber, before use, ensure that the granules inside the absorber are settled and level, this is to provide optimum performance. This can be achieved by a few gentle taps of the absorber between each filling stage.

2. Always refer to the instructions for Use of the anaesthetic workstation concerned.

3. Once opened, any unused content of jericans should be left in sealed jerican and used within one month.

4. Any unused contents of the 1 kg bag should be discarded.

5. Once a re-fillable absorber is filled, the contents should be discarded after one month, regardless of degree of usage.

6. When exhausted, the absorbent must be replaced to ensure that CO₂ absorption in the breathing system is continued.

WHEN TO REPLACE THE ABSORBENT

WARNING A new filter or HMEF should be used on the patient Y piece of the breathing system for every patient. When this is done, the absorbent is not restricted to single patient use and can be used until exhaustion is identified. It is advised that a filter or HMEF is also placed at the expiratory port of the machine.

CAUTIONS

1. If the absorbent has been within the absorber for one month, it is recommended to replace it regardless of extent of usage.

2. The degree of colour change at exhaustion depends on the conditions of use and is an indication only. The absorbent must be used in conjunction with inspiratory FiCO₂ monitoring. The absorbent needs replacing when the inspiratory FiCO₂ concentration exceeds 5mmHg (approximately 0.5% volume FiCO₂).

NOTE: With some absorbents, the colour indicator can sometimes regenerate slowly in between uses and the absorbent may revert back to its original colour. However, this is not a regeneration of the absorption capacity. Therefore, the Intersurgical absorbent should be discarded in the first instance that exhaustion has been determined. Always refer to and follow the Instructions For Use of the anaesthetic machine being used.

COLOUR CHANGES ARE:

White to Violet for 2173000, 2175000 or 2179000

Pink to White for 2172000, 2174000 or 2180000

Green to Violet for 2178000

SAFETY INFORMATION

The absorbent may be used with oxygen, nitrous oxide, Xenon, Enflurane, Halothane, Isoflurane, Desflurane and Sevoflurane.

CAUTIONS

1. Do not use the absorbent if the packaging is breached or damaged, or if there is any evidence of tampering, or if there is any abnormal appearance of the absorbent granules.
2. Before using the absorbent it is essential to become familiar with and follow the Instructions For Use of the anaesthetic machine on which it is to be used.
3. The absorbent must be used in conjunction with monitoring of FiCO_2 inspiratory carbon dioxide.
4. Once exhausted, the absorbent cannot be reused. Do not leave within the machine once exhaustion is identified.

WARNINGS

1. Do not use with trichloroethylene or chloroform because toxic products may be produced by chemical reactions.
2. The absorbent has a moisture content set in the manufacturing process. Do not flush the absorbent with dry gas for long periods. This may result in loss of moisture and reduced CO_2 absorption efficiency. Excessive loss of moisture can lead to failure to absorb CO_2 and to clinical problems.
3. Do not sterilize, disinfect or clean the absorbent before or during use.
4. Do not use the absorbent after expiry date.
5. Never allow pure carbon dioxide to pass through any Carbon Dioxide Absorbents, as this will generate high temperatures.
6. Do not use the absorbent if it becomes excessively hot.
7. Do not use the absorbent if it is believed that dry gas has been left running during long periods of non-use at high flows.
8. Do not use if there is a delay in the increase of inspired agent concentration.
9. It is recommended to use water traps in the patient circuit during low flow anaesthesia (fresh gas flows of 1 L/min or less) and when excessive condensed moisture accumulates in the limbs of the breathing system. Use water traps in both the expiratory and inspiratory limbs to remove water condensation.

For more information please refer to the Intersurgical Absorbents Catalogue or contact Intersurgical.

RISK AND SAFETY PHRASES RELATING TO THE ABSORBENTS WITHIN THE ABSORBENT

H315: Causes skin irritation.

H318 Causes serious eye damage.

P280 Wear protective gloves/protective clothing/eye protection/face protection.

P302/P352 IF ON SKIN: Wash with plenty of soap and water.

P305/351/338. IF IN EYES: Rinse cautiously with water for several minutes.

Remove contact lenses, if present and easy to do. Continue rinsing.

P332/313: If skin irritation occurs: Get medical advice/attention.

STORAGE

The absorbent must be stored, sealed, in a clean, dry environment away from direct sunlight and not under heavy load.

Recommended storage at room temperature over the indicated product shelf-life.

DISPOSAL

Intersurgical advises that after clinical use, the waste absorbent is treated at source with minimal handling; by incineration within a mixture of other clinical waste. In any case, the absorbent must be disposed of according to local regulations. Refer to the MSDS for further details.

DATA SHEET

The Material Safety Data Sheet is available on request.

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Intersurgical Ireland

Tel. +353 (1) 6978540
e-mail: info@intersurgical.ie
249 Blanchardstown Corporate Park 2,
Blanchardstown, Dublin 15, Ireland

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Rx ONLY

INTERSURGICAL®
COMPLETE RESPIRATORY SYSTEMS

Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 9656 300 info@intersurgical.com www.intersurgical.com

Ireland
T: +353 (0)1 697 8540
info@intersurgical.ie

France
T: +33 (0)1 48 76 72 30
info@intersurgical.fr

Deutschland
T: +49 (0)2241 25690
info@intersurgical.de

España
T: +34 91 665 73 15
info@intersurgical-es.com

Portugal
T: +351 219 108 550
info@intersurgical.pt

Italia
T: +39 0535 20836
info@intersurgical.it

Nederland
T: +31 (0)413 243 860
info@intersurgical.nl

Belgie/Belgique
T: +32 5 2456991
info@intersurgical.be

Luxembourg
T: +32 5 2456991
info@intersurgical.be

Sverige
T: +46 08 514 30 600
info@intersurgical.se

Danmark
T: +45 28 69 80 00
info@intersurgical.dk

Lietuva
T: +370 387 66611
info@intersurgical.lt

Россия
T: +7 495 7716809
info@intersurgical.ru

Česká Republika
T: +420 272 911 814
info@intersurgical.cz

Türkiye
T: +90 216 468 88 28
info@intersurgical.com.tr

South Africa
T: +27 (0)11 444 7968
info@intersurgical.co.za

中国
T: +86 519 69817000
info@intersurgical-cn.com

日本
T: 81 (0)3 6863 4359
info@intersurgical.co.jp

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Philipinas
T: +63 2820 4124
info@intersurgical.ph

USA
T: +1 800 828 9633
support@intersurgicalinc.com

Canada
T: +1 905-319-6500
info@intersurgical.ca

Colombia
T: +57 1 742 5161
info@intersurgical.com

Australia
T: +61 (0)2 8048 3300
info@intersurgical.com.au

Distributed in the USA by Intersurgical Incorporated, 6757 Kinne Street, East Syracuse, NY 13057

T: (800) 828-9633 F: (315) 451-3696 support@intersurgicalinc.com www.intersurgical.com

Absorbents-3-F 7005 Issue 7 RO 08.21