10mm neonatal heated wire breathing system

en	10mm neonatal heated wire breathing system
fr	Circuit respiratoire 10 mm de néonatalogie avec fil chauffant
de	Beheiztes Beatmungssystem mit 10-mm-Schlauch für Neugeborene
es	Sistema de respiración con cable caliente de 10 mm para uso neonatal
pt	Sistema de respiração de fio aquecido de 10 mm para doentes recém-nascidos
it	Sistema di respirazione con cavo riscaldato neonatale da 10 mm
nl	10 mm neonataal ademhalingssysteem met verwarmde draad
no	10 mm pustesystem med varmetråd for nyfødte
fi	10 mm:n lämmitysjohtimellinen hengitysjärjestelmä vastasyntyneille
sv	10 mm andningssystem med värmeslinga för nyfödda
da	10 mm neonatal respirations-system med varmetråd
el	Σύστημα αναπνοής με θερμαινόμενο σύρμα 10 mm για νεογνά
lt	10 mm naujagimiams skirta kvėpavimo sistema su šildoma viela
pl	System oddychania o średnicy 10 mm z drutem grzewczym dla noworodków
ru	Дыхательная система 10 мм с подогревом для новорожденных
cs	10mm dýchací systém s vyhřívacím drátem pro novorozeneckou péči
hu	10 mm-es fűtőcsőves újszülött lélegeztető rendszer
sl	Neonatalni dihalni sistem z grelno žico, 10 mm
lv	10 mm apsildāmā vada elpošanas sistēma jaundzimušiem pacientiem
et	10 mm vastsündinute hingamissüsteem soojendusjuhtmega
bg	Неонатална система за дишане с нагревател 10 mm
sr	Sistem za disanje od 10 mm sa grejnim navojem za neonatalne pacijente
ro	Circuit ventilație pentru nou-născuți, cu fir încălzit de 10mm
sk	10 mm dýchací systém s ohrievaným drôtom pre novorodencov
hr	Sustav za disanje za novorođenčad sa zagrijanom žicom od 10 mm
tr	10 mm ısıtma telli yenidoğan solunum sistemi
ja	10mm新生児用熱線入り呼吸回路
zh-s	10 mm 新生儿加热丝呼吸系统
zh-t	10 mm 新生兒加熱線呼吸系統
ar	نظام تنفس مزود بسخان ساكني لحديثي الولادة مقاس 10 مم

CAUTION This product is not suitable for use in an MRI environment.

WARNING To ensure compatibility with this breathing system, only use ventilators that are compliant to ISO 80601-2-12 or ISO 80601-2-72. Prior to use consult the respective documentation and instructions for use of all connected devices, or combination of devices and perform system tests as required.

WARNING

1. Disconnect breathing system from patient during defibrillation.
2. The materials used in the airways and related equipment may not be compatible with anaesthetic or respirable gases, solutions/suspensions/emulsions that have not been evaluated.
3. The use of nebulised drugs and humidification may increase the expiratory filter resistance to flow (where used). Frequently monitor resistance and replace as required.
4. Ventilator accuracy can be affected by use of a gas driven nebuliser.
5. Adding attachments or other devices to the breathing system can change the pressure gradient across the breathing system and adversely affect the ventilator performance. Perform a self-test of the ventilator after full assembly of the breathing system and all connected devices before use on each patient.
6. No modification of this equipment is allowed.
7. Do not use the device past the expiry date.
8. Monitor breathing system condensation and clear any mobile condensate that may build up in the tubes during treatment. Ensure condensate is drained towards the machine and away from the patient.

CAUTION

1. This product has been verified for a maximum duration of 7 days use.
2. Interference can occur between humidifier electrical devices in close proximity. Ensure normal operation of devices prior to use.

Patient category

Intended for use with neonatal patients.

Instructions for use

Pre-use Checks

WARNING

1. Every component of the breathing system must be visually inspected and checked for function, leakage and occlusions, immediately before use on each patient. Ensure yellow that connections are secure.
2. To prevent disconnection of the tubing system during use, use only with devices, and/or accessories in compliance with the applicable international standards.
3. Perform a self-test of the ventilator after full assembly of the breathing system and associated accessories before use on each patient.
4. Prior to use, consult the respective documentation and instructions for use of all connected devices, or combination of devices.
5. Check coiled heated wire. Do not use heated wire breathing system if wire is bunched within the tubes.
6. Ensure that the humidifier is mounted at a lower level than the patient and ventilator. See figure 1.
7. This heated wire breathing system is not suitable for use with flows less than 5L/min.

CAUTION

1. When connecting or disconnecting the system, hold connector and use a push and twist action.
2. Use the red dust cap to occlude the breathing system patient connection port as and when required in the ventilator self-test protocol. Ensure red dust cap is held securely to prevent unwanted disconnection during self-test.

Preparation

1. Fit the chamber

- Slide the humidification chamber on to base unit top plate and remove cassette from chamber ports.
- 2. **Hang water bag from pole**
 - Unwind the fill set from the cassette and remove the red cap from spike, then spike the bag. Hang bag on pole above the humidification chamber. See humidification chamber IFU for details
 - Ensure that the feed line is not kinked and chamber begins to fill.

3. Connect the breathing system (see figure 1)

- Attach dry limb (A) to ventilator inspiratory port and to chamber port.
- Connect the elbow on the green inspiratory limb (B) to other port on chamber.
- Connect expiratory limb (C) to expiratory port on ventilator.
- Attach pressure monitoring line (E) to the ventilator pressure monitoring port (where applicable).

4. Connect temperature probe

- Connect the temperature probe connection to the corresponding socket on

humidifier base unit.

- Insert the distal temperature probe into the port on the elbow at the chamber.

See figure 2

- Insert the proximal temperature probe into the port at the patient end of the green inspiratory limb. See figure 3.

WARNING Ensure proximal temperature probe is located in the inspiratory limb and not in the y-piece.

5. Connect heater wire electrical adapter

- Connect the heated wire adapter into the socket on humidifier base unit.
- Connect the clover leaf end to the socket on the breathing system elbow above chamber.

Dual heated wires breathing systems

- Connect the oval end to socket on elbow at ventilator expiratory port.

6. Turn on humidifier base unit

- Ensure correct mode is activated.

7. Connect breathing system

- Remove the yellow dust cap (D) from y-piece and attach y-piece to patient interface.

In-use checks

WARNING

1. Do not cover the breathing system with anything that may insulate the air flow around the tube as this may impede normal heat transfer and cause damage to the system or affect gas delivery to the patient.
2. The 25cm of breathing system from the patient connection is an applied part, avoid extended contact with patient's skin for the remainder of the system.
3. Do not stretch the tubing.
4. Monitor the clinical parameters of the patient during the treatment at regular intervals.

CAUTION

1. Check humidification chamber water level before starting and during treatment. See chamber IFU for more details.
2. Check humidifier status before starting and during treatment. See humidifier IFU for more details.

Duration of use

Single use product. Maximum duration of use 7 days. The start of use is defined from the first connection to the equipment. Do not reuse. Reuse poses a threat to patient safety in terms of cross infection and device malfunction and breakage. Cleaning can lead to device malfunction and breakage.

Note

WARNING Do not use outside the specified environmental conditions. Using the device outside this range can affect the quality of the therapy or injure the patient.

Ensure compatibility of all parts of the system when connecting accessories and breathing systems to a ventilator and/ or humidifier.

Technical and performance data

WARNING Ensure that the breathing system technical and performance data are appropriate to the intended patient and treatment.

Patient ventilation and vital signs should be monitored. Please refer to product label for full information.

Humidity output:

- **Invasive mode:** ≥ 33mg/L* @ 5-15L/min
- **Non-invasive mode:** ≥ 12mg/L* @ 5-15L/min

Working pressure: <125 cmH₂O

Maximum operating temperature: 25 °C

Minimum operating temperature: 18 °C

Breathing system leakage: < 70mL/min, additional limbs < 25mL/min at 60cmH₂O.

* *At sea level when used with Intersurgical humidification chamber.*

Recommended storage at room temperature over the indicated product shelf-life.

For further information contact Intersurgical - info@intersurgical.com.

Disposal

Following use, the product must be disposed of in accordance with local hospital, infection control and waste disposal regulations.

Standards

This product complies with the following relevant standards:

- BS EN ISO 5367: Anaesthetic and respiratory equipment. Breathing tube and sets
- BS EN ISO 5356-1: Anaesthetic and respiratory equipment. Conical connectors. Cones and sockets
- ISO 80601-2-74: Medical electrical equipment - Particular requirements for basic safety and essential performance of respiratory humidifying equipment
- BS EN ISO 80601-2-12:2011 MEE Critical Care Ventilators.

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INTERSURGICAL®
COMPLETE RESPIRATORY SYSTEMS

Intersurgical Ltd, Crane House, Molly Millars Lane,
Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 9656 300 info@intersurgical.com
www.intersurgical.com

Ireland

T: +353 (0)1 697 8540
info@intersurgical.ie

France

T: +33 (0)1 48 76 72 30
info@intersurgical.fr

Deutschland

T: +49 (0)2241 25690
info@intersurgical.de

España

T: +34 91 665 73 15
info@intersurgical-es.com

Portugal

T: +351 219 108 550
info@intersurgical.pt

Italia

T: +39 0535 20836
info@intersurgical.it

Nederland

T: +31 (0)413 243 860
info@intersurgical.nl

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T: +32 5 2456991
info@intersurgical.be

Luxembourg

T: +32 5 2456991
info@intersurgical.be

Sverige

T: +46 08 514 30 600
info@intersurgical.se

Danmark

T: +45 28 69 80 00
info@intersurgical.dk

Lietuva

T: +370 387 66611
info@intersurgical.lt

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T: +7 495 7716809
info@intersurgical.ru

Česká Republika

T: +420 272 911 814
info@intersurgical.cz

Türkiye

T: +90 216 468 88 28
info@intersurgical.com.tr

South Africa

T: +27 (0)11 444 7968
info@intersurgical.co.za

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T: +86 519 69817000
info@intersurgical-cn.com

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T: 81 (0)3 6863 4359
info@intersurgical.co.jp

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T: +886 4 2380 5430
info@intersurgical.com.tw

Pilipinas

T: +63 2820 4124
info@intersurgical.ph

USA

T: +1 800 828 9633
support@intersurgicalinc.com

Canada

T: +1 905-319-6500
info@intersurgical.ca

Colombia

T: +57 1 742 5161
info@intersurgical.com

Australia

T: +61 (0)2 8048 3300
info@intersurgical.com.au

Dual float auto fill humidification chamber

en	DUAL FLOAT AUTO FILL HUMIDIFICATION CHAMBER
fr	CHAMBRE D'HUMIDIFICATION À REMPLISSAGE AUTOMATIQUE AVEC DEUX FLOTTEURS
de	BEFEUCHTERKAMMER MIT AUTOMATISCHER FÜLLUNG UND DOPPELSCHWIMMER
es	CÁMARA DE HUMIDIFICACIÓN DE LLENADO AUTOMÁTICO CON DOBLE FLOTADOR
pt	CÂMARA DE HUMIDIFICAÇÃO DE ENCHIMENTO AUTOMÁTICO COM FLUTUADOR DUPLO
it	CAMERA DI UMIDIFICAZIONE A RIEMPIMENTO AUTOMATICO A DOPPIO GALLEGGIANTE
nl	AUTOMATISCHE VUL BEVOCHTIGINGSKAMER MET DUBBELE VLOTTER
no	FUKTINGSKAMMER MED DOBBEL FLOTTØR OG AUTOMATISK FYLING
fi	AUTOMAATTISESTI TÄYTETTÄVÄ KAKSOISVENTTIILILLÄ VARUSTETTU KOSTUTUSKAMMIO
sv	BEFUKTNINGSKAMMARE MED DUBBELT FLÖDE OCH AUTOMATISK Fyllning
da	BEFUGTNINGSKAMMER MED TO FLYDERE, DER PÅFYLDES AUTOMATISK
el	ΘΑΛΑΜΟΣ ΥΓΡΑΝΣΗΣ ΑΥΤΟΜΑΤΗΣ ΠΛΗΡΩΣΗΣ ΜΕ ΔΙΠΛΟ ΠΛΩΤΗΡΑ
lt	DVIAUKŠTĖ AUTOMATINIO UŽPILDYMO DRĖKINIMO KAMERA
pl	AUTOMATYCZNE NAPEŁNIANA KOMORA NAWILŻANIA Z DWOMA PŁYWAKAMI
ru	САМОЗАПОЛНЯЮЩАЯСЯ КАМЕРА УВЛАЖНЕНИЯ С ДВОЙНЫМ ПОПЛАВКОМ
cs	AUTOMATICKY PLNĚNÁ ZVLHČOVACÍ KOMORA SE SYSTÉMEM DVOU PLOVÁKŮ
hu	KETTŐS, ÚSZÓVAL ELLÁTOTT AUTOMATIKUS FELTÖLTÉSŰ PÁRÁSÍTÓ KAMRA
sl	VLAŽILNA KOMORA Z DVOJNIM PLOVCEM IN SAMODEJNIM POLNJENJEM
lv	AUTOMĀTISKI UZPILDĀMA MITRINĀŠANAS KAMERA AR DIVIEM PLUDIŅIEM
et	KANE UJUKIGA AUTOMAATSE TÄITUMISEGA NIISUTIKAMBER
bg	ОВЛАЖНИТЕЛНА КАМЕРА С АВТОМАТИЧНО ПЪЛНЕНЕ С ДВОЕН ПОПЛАВЪК
sr	KOMORA ZA OVLAŽIVANJE SA DVA VENTILA I AUTOMATSKIM PUNJENJEM
ro	CAMERĂ DE UMIDIFICARE CU UMPLERE AUTOMATĂ CU FLOTOARE DUALE
sk	ZVLHČOVACIA KOMORA S DVOJPLAVÁKOVÝM SYSTÉMOM AUTOMATICKÉHO PLNENIA
hr	KOMORA ZA OVLAŽIVANJE S AUTOMATSKIM PUNJENJEM I DVA PLOVKA
tr	ÇİFT ŞAMANDIRALI, OTOMATİK DOLUMLU NEMLENDİRME HAZNESİ
ja	デュアルフロート自動充填加湿器チャンバ
zh-s	双浮子自动加水式湿化罐
zh-t	雙浮式自動填充加湿瓶
ar	غرفة ترطيب تلقائية الملء مزدوجة العوامة

The instruction for use is not supplied individually. Only 1 copy of the instructions for use is provided per box and should therefore be kept in an accessible location for all users. More copies are available upon request.

NOTE: Distribute instructions to all product locations and users.

These instructions contain important information for the safe use of the product. Read these instructions for use in their entirety, including warnings and cautions before using this product. Failure to properly follow warnings, cautions and instructions could result in serious injury or death to the patient. Serious incidents to be reported as per the local regulation.

WARNING: A WARNING statement provides important information about a potentially hazardous situation which, if not avoided could result in death or serious injury.

CAUTION: A CAUTION statement provides important information about a potentially hazardous situation which, if not avoided could result in minor or moderate injury.

Description. Auto fill humidification chamber with fill set. **Intended use.** To humidify the respiratory gases being delivered to a patient.

Intended user. For use by qualified medical personnel. This product may be used in a homecare setting. Appropriately qualified medical personnel to instruct homecare users on the correct usage.

Indications

For use in a hospital environment.

For use in a homecare environment.

For use with the Intersurgical PMH7000, Fisher & Paykel® MR850 and MR730 humidifiers.

* Fisher and Paykel is a registered trademark of Fisher & Paykel® Healthcare Limited.

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Patient category. Intended for use with adult, paediatric and neonatal patients.

INSTRUCTIONS FOR USE

Pre-use Checks

CAUTION When connecting or disconnecting the system, use a push and twist action.

WARNING To prevent disconnection of the tubing system during use, use only tubes in compliance with applicable ISO standards.

WARNING Every component of the breathing system must be visually inspected and checked for function, leakage and occlusions, immediately before use on each patient. Ensure that connections are secure.

WARNING Perform a self-test of the ventilator after full assembly of the breathing system and associated accessories before use on each patient.

WARNING Prior to use, consult the respective documentation and instructions for use of all connected devices, or combination of devices.

WARNING Ensure humidifier is level, maximum tilt angle < 10 degrees. See figure a.

WARNING Only use the chamber with sterile water suitable for inhalation or equivalent.

Preparation

1. Fit the chamber

• Slide the humidification chamber on to base unit top plate and remove the green cassette from the humidification ports. See figure a - 3.

2. Hang water bag from pole

• Unwind the fill set from the cassette and remove the red cap from spike (see figure a - 1), then spike the bag (see figure a - 2). Hang bag so that the height (H)

above the chamber is greater than maximum ventilated pressure (Pmax) in cmH₂O. See figure a.

• Ensure that the feed line is not kinked and chamber begins to fill.

• Open spike vent when using rigid water reservoir. See figure a - 4.

CAUTION: Ensure that the feed line is free from trapped air bubbles before commencing treatment.

3. Turn on humidifier base unit

• Ensure correct mode is activated.

In-use checks

CAUTION Check humidification chamber water level before starting and during treatment. See figure a.

CAUTION Check humidifier status before starting and during treatment. See humidifier IFU for more details.

WARNING Replace the humidifier if water level exceeds maximum fill line. See figure a.

CAUTION When using Infant nasal CPAP please use a pressure infuser to ensure that reservoir bag pressure exceeds 140 cmH₂O. Check fluid levels after 2 minutes and while in use.

WARNING Metal surfaces may be hot when exposed. See figure a.

WARNING Monitor the clinical parameters of the patient, i.e O₂ concentration measured at the point of delivery to the patient, during the treatment.

Duration of use

Single use product. Maximum duration of use 7 days. Do not reuse. Reuse poses a threat to patient safety in terms of cross infection and device malfunction and breakage. Cleaning can lead to device malfunction and breakage.

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Environmental Operating Conditions: Performance verified at room/ ambient temperature at sea level.

Resistance to flow

	Resistance to Flow (cm H ₂ O)					
	2.5 L/min	15 L/min	30 L/min	60 L/min	90 L/min	120 L/min
Max Fill Chamber	0.0	0.0	0.0	0.4	0.9	1.5

Compliance

	Compliance at 60cmH ₂ O
Min Fill Chamber	0.5
Max Fill Chamber	0.3

Min Humidification output: 33 mg/L @ BTPS

WARNING Please refer to breathing system IFU for humidification output data. Ensure technical and performance data are appropriate to the intended patient and treatment.

Patient ventilation and vital signs should be monitored.

For further information contact Intersurgical – info@intersurgical.com.

CAUTION Only to be used by qualified medical personnel or a user suitably trained to use the device by qualified medical personnel.

CAUTION Do not stretch or cut the fill set.

CAUTION This product has been verified for a maximum duration of 7 days use.

CAUTION Following use, the product must be disposed of in line with the local healthcare, hygiene and waste disposal protocols and regulations.

WARNING The materials used in the airways and related equipment may not be compatible with anaesthetic or respirable gases, solutions/suspensions/emulsions that have not been evaluated.

WARNING The use of nebulised drugs may cause discolouration of the chamber and / or water.

WARNING Do not use sterile products if the packaging is damaged, refer to the product label for information.

WARNING Use of this breathing system with systems and accessories not compliant with applicable ISO standards can result in degraded performance. Prior to use, consult the respective documentation and instructions for use of all connected devices, or combination of devices to ensure

compatibility.

WARNING Do not use outside the specified environmental conditions, as per the Humidifier IFU. Using the device outside this range can affect the quality of the therapy or injure the patient.

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- BS EN ISO 5356-1: Anaesthetic and respiratory equipment. Conical connectors. Cones and sockets.
- ISO 80601-2-74: Medical electrical equipment – Particular requirements for basic safety and essential performance of respiratory humidifying equipment.
- BS EN ISO 18190: Anaesthetic and respiratory equipment. General requirements for airways and related devices.

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INTERSURGICAL[®]
COMPLETE RESPIRATORY SYSTEMS

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Česká Republika
T: +420 272 911 814
info@intersurgical.cz

Türkiye
T: +90 216 468 88 28
info@intersurgical.com.tr

South Africa
T: +27 (0)11 444 7968
info@intersurgical.co.za

中国
T: +86 519 69817000
info@intersurgical-cn.com

日本
T: 81 (0)3 6863 4359
info@intersurgical.co.jp

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Pilipinas
T: +63 2820 4124
info@intersurgical.ph

USA
T: +1 800 828 9633
support@intersurgicalinc.com

Canada
T: +1 905-319-6500
info@intersurgical.ca

Colombia
T: +57 1 742 5161
info@intersurgical.com

Australia
T: +61 (0)2 8048 3300
info@intersurgical.com.au

patent: <http://us.intersurgical.com/support>

MADE BY INTERSURGICAL IN LITHUANIA

Distributed in the USA by Intersurgical Incorporated, 6757 Kinne St., East Syracuse, NY 13057
T: (800) 828-9633 F: (315) 451-3696 support@intersurgicalinc.com www.intersurgical.com

Neonatal water trap

en	Neonatal water trap
fr	Piège à eau pour nouveau-nés
de	Wasserauffangbehälter für Neugeborene
es	Colector de agua neonatal
pt	Separador de água neonatal
it	Separatore di condensa per uso neonatale
nl	Neonatale watervanger
no	Neonatal vannfelle
fi	Vastasyntyneille tarkoitettu vedenerotin
sv	Neonatal vattenavskiljare
da	Neonatal vandfælde
el	Υδατοπαγίδα νεογνών
lt	Naujagimių vandens skirtuvas
pl	Skraplacz dla noworodków
ru	Неонатальный водоотделитель
cs	Neonatální kondenzační nádobka
hu	Újszülöttekhez szolgáló vízleválasztó
sl	Neonatalni vodni zbiralnik
lv	Neonatālais ūdens uztvērējs
et	Veeseparaator kasutamiseks vastsündinutel
bg	Неонатален водоуловител
sr	Uređaj za sakupljanje vode tokom disanja na respiratoru kod novorođenčad
ro	Separator de apă neonatal
sk	Neonatálny odlučovač vody
hr	Neonatalni kolektor vode
tr	Yenidoğan su tutucusu

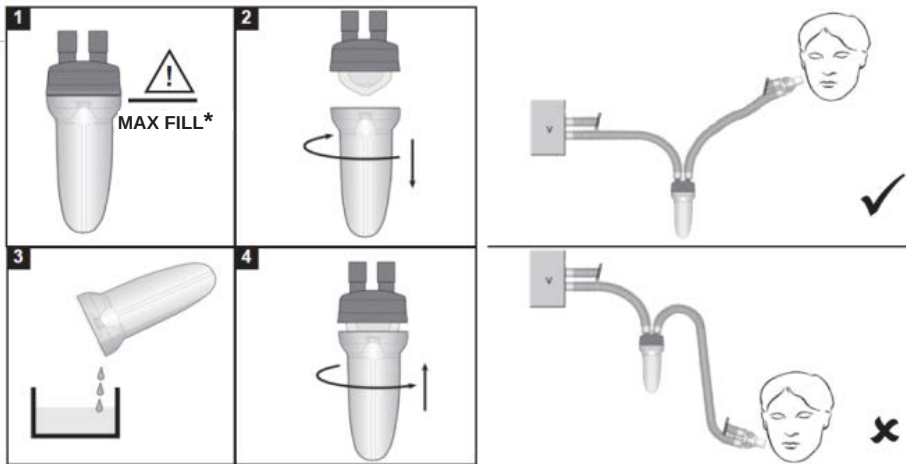
Figure 1. Correct set-up of water trap/ Figure 1. Installation correcte du piège à eau/ Abbildung 1. Korrekte Einrichtung des Wasserauffangbehälters/ Figura 1. Preparación correcta del colector/ Figura 1. Configuração correta do separador de água/ Figura 1. Posizionamento corretto del separatore di condensa/ Afbeelding 1. Juiste opstelling van de watervanger/ Figur 1. Riktig oppsett av vannfellen/ Kuva 1. Vedenerottimen oikea asennus/ Figur 1. Korrekt montering av vattenavskiljaren/ Figur 1. Korrekt opsætning af vandfælde/ Εικόνα 1. Σωστή εγκατάσταση της υδατοπαγίδας/ 1 pav. Tinkamas vandens skirtuvo nustatymas/ Rysunek 1. Poprawna konfiguracja skraplacza/ Рис. 1. Правильная установка водоотделителя/ Obrázek 1. Správná instalace kondenzační nádoby/ 1. ábra: A vízelválasztó helyes beállítása/ Slika 1. Pravilna namestitev vodnega zbiralnika/ 1. attēls. Pareiza ūdens uztvērēja uzstādīšana/ Joonis 1. Veeseparaatori korrektn seadistamine/ Фигура 1. Правилна постановка на водоуловителя/ Slika 1. Ispravno postavljanje uređaja za sakupljanje vode/ Figura 1. Instalarea corectă a separatorului de apă/ Obrázok č. 1. Správne nastavenie odlučovača vody/ Slika 1. Pravilno postavljanje kolektora vode/ Şekil 1. Su tutucusunun doğru kurulumu

* REMPLISSAGE MAXIMUM/ MAX FÜLLSTAND/ RELLENO MÁX/ ENCHIMENTO MÁXIMO/ LIVELLO MAX/ MAX. VULLING/ MAKS. FÜLLING/ MAKSIMITÄYTTÖ/ MAXIMAL FÜLLNAD/ MAKS. FYLDNING/ ΜΕΓ. ΠΛΗΡΟΣΗ/ MAKS. UZPILDYMAS/ MAKS. POZIOM NAPELNIENIA/ MAKS. ЗАПОЛНЕНИЕ/ MAX. HLADINA/ MAX. SZINT/ MAKS. NAPOLNENOST/ MAKSIMĀLĀ AIZPILDE/ MAKSIMAALNE TĀITUVUS/ МАКСИМАЛНО НАПЪЛВАНЕ/ MAKS. PUNJENJE/ UMPLERE MAXIMĂ/ MAXIMUM/ MAKS. PUNJENJE/ MAKSIMUM DOLUM

en

The instruction for use (IFU) is not supplied individually. Only 1 copy of the instructions for use is provided per box and should therefore be kept in an accessible location for all users. More copies are available upon request.

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Description. Neonatal water trap.

Intended use.

To collect condensed water in a breathing system.

Intended user

For use by qualified medical personnel. This product may be used in a homecare setting. Appropriately qualified medical personnel to instruct homecare users on correct usage.

Indications

This water trap is intended to be used with Intersurgical breathing systems.

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WARNING Perform a self-test of the ventilator after full assembly of the breathing system and associated accessories

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WARNING Ensure water trap is positioned at the lowest point and is level. See figure 1.

CAUTION When connecting or disconnecting the water trap reservoir, use a locate and twist action. See figure 1.

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In-use checks

WARNING Ensure excess fluid in the system is collected in the water trap reservoir.

WARNING Monitor fluid level in water trap reservoir during treatment and empty before completely full.

CAUTION Empty trap contents in line with hospital contaminated liquid waste protocol.

WARNING Monitor the clinical parameters of the patient during the treatment.

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- BS EN ISO 5367: Anaesthetic and respiratory equipment. Breathing tube and sets.
- BS EN ISO 5356-1: Anaesthetic and respiratory equipment. Conical connectors. Cones and sockets.
- BS EN ISO 18190: Anaesthetic and respiratory equipment. General requirements for airways and related devices.

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