

L A N D A N G E R	FICHE TECHNIQUE PRODUIT PRODUCT TECHNICAL SHEET	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-------------------	--	--

Editée le / Publishing date

2022-04-28

Page

1

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
C010001	CROCHET A PHLEBECTOMIE DU DR MULLER PTE FINE RECOURBEE DROITIER N° 1 L:150MM MODELE LEPINE	3661219291800
	DR MULLER PHLEBECTOMY HOOK - FINE POINT CURVED RIGHT-HANDED N° 1 L:150MM - LEPINE MODEL	 3 661219 291800

Catalogues / catalogs

Landanger 2015

Chapitre-page

05-007

DESCRIPTIF DU DISPOSITIF/ *PRODUCT DESCRIPTION*

L A N D A N G E R

Phlébectomie - Phlebectomy - Flebectomía

TULIPE



1.00
C012000
12 cm - 4 3/4"

1.25
C012001
12 cm - 4 3/4"

1.50
C012002
12 cm - 4 3/4"


1.75
C012003
12 cm - 4 3/4"

2.00
C012004
13 cm - 5"

MULLER

Crochets à phlébectomie - Phlebectomy hooks - Ganchos para flebectomía

N° 1	C010106	15 cm - 6"	C010001	N° 1
N° 2	C010206	15 cm - 6"	C010002	N° 2
N° 3	C010306	15 cm - 6"	C010003	N° 3
N° 4	C010406	15 cm - 6"		

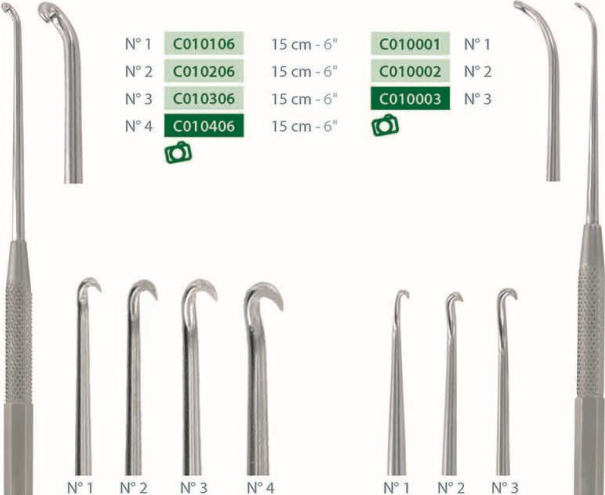

N° 1 N° 2 N° 3 N° 4 N° 1 N° 2 N° 3

Modèle Lepine
Lepine model
Modelo Lepine



RUMEL BELMONT

TOURNIQUET VASCULAIRE
VASCULAR TOURNIQUET
TOURNIQUETE VASCULAR

C09150
32 cm - 12 1/2"

Mandrin seul
Single pad
Mandril único

C091510
26 cm - 11"

V-Clamps
V-Clamps
V-Clamps

V-7

Copyright © LANDANGER - 2009 - 2015

LANDANGER	FICHE TECHNIQUE PRODUIT <i>PRODUCT TECHNICAL SHEET</i>	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-----------	---	---

Éditée le / *Publishing date*

2022-04-28

Page

2

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
C010001	CROCHET A PHLEBECTOMIE DU DR MULLER PTE FINE RECOURBEE DROITIER N° 1 L:150MM MODELE LEPINE	3661219291800
	DR MULLER PHLEBECTOMY HOOK - FINE POINT CURVED RIGHT- HANDED N° 1 L:150MM - LEPINE MODEL	

INFORMATIONS REGLEMENTAIRES / *REGULATORY INFORMATION*

Classe CE :

I

Code GMDN :

46768

CE Class :

GMDN code :

PERFORMANCES REVENDIQUEES / *INTENDED PERFORMANCE*

Dispositif médical de bloc opératoire :

REUTILISABLE/REUSABLE

Operating room medical device :

CARACTERISTIQUES SPECIFIQUE / *SPECIFIC FEATURES*

Type matière :

ACIER INOXYDABLE / STAINLESS STEEL

Material type :

Nuance matière :

X30Cr13(1.4028) + X2CrNi19-11(1.4306)

Material grade :

Dimensions :

Voir désignations produits

Landanger 2015

Outer dimensions :

See products description

Chapitre-pages

05-007

Nettoyage/Pré-désinfection : Standard pour instrumentation de bloc opératoire réutilisable.

Cleaning / Pre-disinfection : *Standard conditions for reusable surgical instruments*

Stérilisation :

La stérilisation doit être effectuée selon un procédé de stérilisation à la vapeur conformément à la norme ISO 17665-1 (procédé sous vide fractionné) avec un stérilisateur conformément à la norme EN 285.

Méthode de stérilisation à la vapeur saturée	Température	Durée	Temps de séchage	Pression absolue
Pré-vide	132 °C (270 °F)	4 minutes	30 minutes	2,9 bars, 2,8 atm, 41,6 psi
Pré-vide	134 °C (273 °F)	3 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	5 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	18 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi

Sterilization process :

The Sterilization shall be performed according to ISO 17665-1 validated steam sterilisation procedure (fractional vacuum procedure) in an autoclave in compliance with EN 285.

Saturated steam sterilization method	Temperature	Time	Drying time	Absolute pressure
Pre-vacuum	132 °C (270 °F)	4 minutes	30 minutes	2.9 bar, 2.8 atm, 41.6 psi
Pre-vacuum	134 °C (273 °F)	3 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	5 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	18 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi

EMBALLAGE / *PACKAGING*

Conditionnement :

Pochette plastique en polyéthylène (PE).

Packaging :

Polyethylene plastic pouch (PE).