

L A N D A N G E R	FICHE TECHNIQUE PRODUIT PRODUCT TECHNICAL SHEET	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-------------------	--	--

Editée le / Publishing date

2022-04-28

Page

1

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
DC70403-25	CANULE D:2,5MM ANGULE 30° A 25MM DU BOUT	3661219319306
	INJECTION CANNULA DIA:2.5MM 30° ANGULATED AT 25MM FROM TIP	 3 661219 319306

Catalogues / catalogs

Delacroix 2020

Chapitre-page

06-006

DESCRIPTIF DU DISPOSITIF/ *PRODUCT DESCRIPTION*




Suction Tubes / Canules d'Aspiration / Canulas de Aspiración



 **DC70401-15**
ø : 1,6 mm
 **DC70401-25**
ø : 2,5 mm

 **DC70402-15**
ø : 1,6 mm
 **DC70402-25**
ø : 2,5 mm

 **DC70403-15**
ø : 1,6 mm
 **DC70403-25**
ø : 2,5 mm

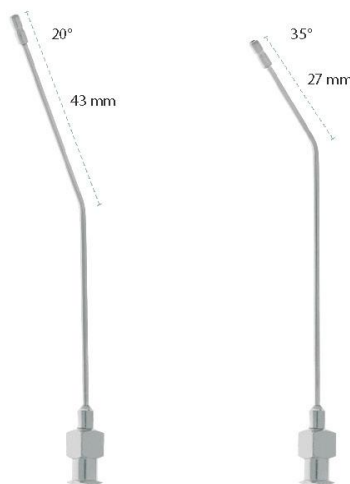



6 - Special

6

Copyright LANDANGER - 2009 - 2020

 1/1

LANDANGER	FICHE TECHNIQUE PRODUIT <i>PRODUCT TECHNICAL SHEET</i>	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-----------	---	---

Éditée le / *Publishing date*

2022-04-28

Page

2

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
DC70403-25	CANULE D:2,5MM ANGULE 30° A 25MM DU BOUT	3661219319306
	INJECTION CANNULA DIA:2.5MM 30° ANGULATED AT 25MM FROM TIP	 3 661219 319306

INFORMATIONS REGLEMENTAIRES / *REGULATORY INFORMATION*

Classe CE :
CE Class :

Ila

Code GMDN :
GMDN code :

57822

PERFORMANCES REVENDIQUEES / *INTENDED PERFORMANCE*

Dispositif médical de bloc opératoire :
Operating room medical device :

REUTILISABLE/REUSABLE

CARACTERISTIQUES SPECIFIQUE / *SPECIFIC FEATURES*

Type matière :
Material type :

ACIER INOXYDABLE / STAINLESS STEEL

Nuance matière :
Material grade :

X8CrNiS18-9(1.4305) + X5CrNi18-10(1.4310) + X2CrNiMo18-15-3(1.4441)

Dimensions :
Outer dimensions :

Voir désignations produits
See products description
Chapitre-pages

Delacroix 2020

06-006

Nettoyage/Pré-désinfection :
Cleaning / Pre-disinfection :

Standard pour instrumentation de bloc opératoire réutilisable.
Standard conditions for reusable surgical instruments

Stérilisation :

La stérilisation doit être effectuée selon un procédé de stérilisation à la vapeur conformément à la norme ISO 17665-1 (procédé sous vide fractionné) avec un stérilisateur conformément à la norme EN 285.

Méthode de stérilisation à la vapeur saturée	Température	Durée	Temps de séchage	Pression absolue
Pré-vide	132 °C (270 °F)	4 minutes	30 minutes	2,9 bars, 2,8 atm, 41,6 psi
Pré-vide	134 °C (273 °F)	3 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	5 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	18 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi

Sterilization process :

The Sterilization shall be performed according to ISO 17665-1 validated steam sterilisation procedure (fractional vacuum procedure) in an autoclave in compliance with EN 285.

Saturated steam sterilization method	Temperature	Time	Drying time	Absolute pressure
Pre-vacuum	132 °C (270 °F)	4 minutes	30 minutes	2.9 bar, 2.8 atm, 41.6 psi
Pre-vacuum	134 °C (273 °F)	3 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	5 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	18 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi

EMBALLAGE / *PACKAGING*

Conditionnement :
Packaging :

Pochette plastique en polyéthylène (PE).
Polyethylene plastic pouch (PE).