

LANDANGER	FICHE TECHNIQUE PRODUIT PRODUCT TECHNICAL SHEET	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-----------	--	---

Éditée le / Publishing date

2022-04-28

Page

1

REFERENCE	DESIGNATIONS - DESCRIPTIONS	CODE EAN
DC82142-30	PAIRE DE VALVES DUBOST REVERSIBLES POUR THORACCESS STERNAL I:50MM P:30MM	3661219292654
	PAIR OF DUBOST UP&DOWN BLADES FOR STERNAL THORACCESS W:50MM D:30MM	3 661219 292654

Catalogues / catalogs

MIS 2021


Chapitre-page

55

DESCRIPTIF DU DISPOSITIF/ PRODUCT DESCRIPTION



SternAccess Retractor / Écarteur SternAccess / Separador SternAccess

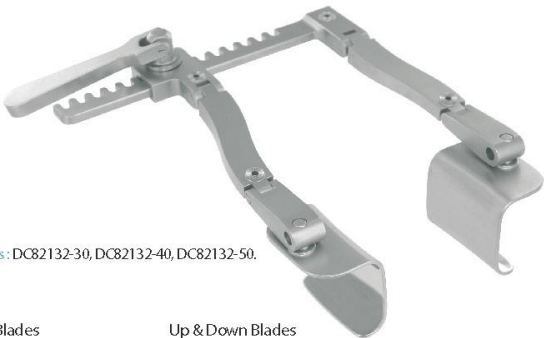

DC82130-00
c : 115 mm - d : 130 mm
Included - Inclus - Incluidos: DC82132-30, DC82132-40, DC82132-50.

Swivel Standard Blades
Valves Standard Rotatives
Valvas Rotativas Estándar




DC82131-30	a : 36 mm - b : 30 mm
DC82131-40	a : 36 mm - b : 40 mm
DC82131-50	a : 36 mm - b : 50 mm
DC82132-30	a : 50 mm - b : 30 mm
DC82132-40	a : 50 mm - b : 40 mm
DC82132-50	a : 50 mm - b : 50 mm
DC82133-30	a : 64 mm - b : 30 mm
DC82133-40	a : 64 mm - b : 40 mm
DC82133-50	a : 64 mm - b : 50 mm

Up & Down Blades
Valves Réversibles
Valvas Reversibles



DC82141-30	a : 36 mm - b : 30 mm
DC82141-40	a : 36 mm - b : 40 mm
DC82141-50	a : 36 mm - b : 50 mm
DC82142-30	a : 50 mm - b : 30 mm
DC82142-40	a : 50 mm - b : 40 mm
DC82142-50	a : 50 mm - b : 50 mm
DC82143-30	a : 64 mm - b : 30 mm
DC82143-40	a : 64 mm - b : 40 mm
DC82143-50	a : 64 mm - b : 50 mm

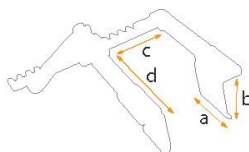

Copyright LANDANGER - 2009 - 2021



Retractors
Écarteurs
Retradores

LANDANGER	FICHE TECHNIQUE PRODUIT <i>PRODUCT TECHNICAL SHEET</i>	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-----------	---	---

Éditée le / *Publishing date*

2022-04-28

Page

2

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
DC82142-30	PAIRE DE VALVES DUBOST REVERSIBLES POUR THORACCESS STERNAL 1:50MM P:30MM	3661219292654
	PAIR OF DUBOST UP&DOWN BLADES FOR STERNAL THORACCESS W:50MM D:30MM	 3 661219 292654

INFORMATIONS REGLEMENTAIRES / *REGULATORY INFORMATION*

Classe CE :

I

Code GMDN :

46774

CE Class :

GMDN code :

PERFORMANCES REVENDIQUEES / *INTENDED PERFORMANCE*

Dispositif médical de bloc opératoire :

REUTILISABLE/REUSABLE

Operating room medical device :

CARACTERISTIQUES SPECIFIQUE / *SPECIFIC FEATURES*

Type matière :

ACIER INOXYDABLE / STAINLESS STEEL

Material type :

Nuance matière :

X30Cr13(1.4028) + X2CrNi19-11(1.4306)

Material grade :

Dimensions :

Voir désignations produits

MIS 2021

Outer dimensions :

See products description

Chapitre-pages

55

Nettoyage/Pré-désinfection : Standard pour instrumentation de bloc opératoire réutilisable.

Cleaning / Pre-disinfection : *Standard conditions for reusable surgical instruments*

Stérilisation :

La stérilisation doit être effectuée selon un procédé de stérilisation à la vapeur conformément à la norme ISO 17665-1 (procédé sous vide fractionné) avec un stérilisateur conformément à la norme EN 285.

Méthode de stérilisation à la vapeur saturée	Température	Durée	Temps de séchage	Pression absolue
Pré-vide	132 °C (270 °F)	4 minutes	30 minutes	2,9 bars, 2,8 atm, 41,6 psi
Pré-vide	134 °C (273 °F)	3 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	5 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	18 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi

Sterilization process :

The Sterilization shall be performed according to ISO 17665-1 validated steam sterilisation procedure (fractional vacuum procedure) in an autoclave in compliance with EN 285.

Saturated steam sterilization method	Temperature	Time	Drying time	Absolute pressure
Pre-vacuum	132 °C (270 °F)	4 minutes	30 minutes	2.9 bar, 2.8 atm, 41.6 psi
Pre-vacuum	134 °C (273 °F)	3 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	5 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	18 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi

EMBALLAGE / *PACKAGING*

Conditionnement :

Pochette plastique en polyéthylène (PE).

Packaging :

Polyethylene plastic pouch (PE).