

L A N D A N G E R	FICHE TECHNIQUE PRODUIT <i>PRODUCT TECHNICAL SHEET</i>	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-------------------	---	---

Éditée le / *Publishing date*

2022-04-28

Page

1

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
J320030	CISEAU OSTEOTOME L:135MM 10MM	3661219347293
	OSTEOTOME L:135MM 10MM	 3 661219 347293

Catalogues / *catalogs*

Landanger 2015

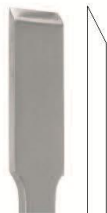
Chapitre-page

13-069

DESCRIPTIF DU DISPOSITIF/ *PRODUCT DESCRIPTION*

L A N D A N G E R

Ciseaux burins, ostéotomes, gouges - Bone chisels, osteotomes, gouges - Osea cincelas, osteotomos, gubias

	Burin Bone chisel Osea cincela	Osteotome Osteotome Osteotomo	Gouge Gouge Gubia
	16 mm J320016	J320036	J360916
	14 mm J320014	J320034	J360914
	12 mm J320012	J320032	J360912
	10 mm J320010	J320030	J360910
	8 mm J320008	J320028	J360908
	6 mm J320006	J320026	J360906
	4 mm J320004	J320024	J360904



J320012
13,5 cm - 5 1/4"





J360912
13,5 cm - 5 1/4"

XIII - 69

XIII - Oseaux
XIII - Bone surgery
XIII - Cirugía osea

Copyright © LANDANGER - 2009 - 2015

LANDANGER	FICHE TECHNIQUE PRODUIT <i>PRODUCT TECHNICAL SHEET</i>	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-----------	---	---

Éditée le / *Publishing date*

2022-04-28

Page

2

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
J320030	CISEAU OSTEOTOME L:135MM 10MM	3661219347293
	OSTEOTOME L:135MM 10MM	 3 661219 347293

INFORMATIONS REGLEMENTAIRES / *REGULATORY INFORMATION*

Classe CE :
CE Class :

I

Code GMDN :
GMDN code :

12844

PERFORMANCES REVENDIQUEES / *INTENDED PERFORMANCE*

Dispositif médical de bloc opératoire :
Operating room medical device :

REUTILISABLE/REUSABLE

CARACTERISTIQUES SPECIFIQUE / *SPECIFIC FEATURES*

Type matière :
Material type :

ACIER INOXYDABLE / STAINLESS STEEL

Nuance matière :
Material grade :

X20Cr13(1.4021) + X5CrNi18-10(1.4301)

Dimensions :
Outer dimensions :

Voir désignations produits
See products description
Chapitre-pages

Landanger 2015

13-069

Nettoyage/Pré-désinfection :
Cleaning / Pre-disinfection :

Standard pour instrumentation de bloc opératoire réutilisable.
Standard conditions for reusable surgical instruments

Stérilisation :

La stérilisation doit être effectuée selon un procédé de stérilisation à la vapeur conformément à la norme ISO 17665-1 (procédé sous vide fractionné) avec un stérilisateur conformément à la norme EN 285.

Méthode de stérilisation à la vapeur saturée	Température	Durée	Temps de séchage	Pression absolue
Pré-vide	132 °C (270 °F)	4 minutes	30 minutes	2,9 bars, 2,8 atm, 41,6 psi
Pré-vide	134 °C (273 °F)	3 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	5 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	18 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi

Sterilization process :

The Sterilization shall be performed according to ISO 17665-1 validated steam sterilisation procedure (fractional vacuum procedure) in an autoclave in compliance with EN 285.

Saturated steam sterilization method	Temperature	Time	Drying time	Absolute pressure
Pre-vacuum	132 °C (270 °F)	4 minutes	30 minutes	2.9 bar, 2.8 atm, 41.6 psi
Pre-vacuum	134 °C (273 °F)	3 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	5 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	18 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi

EMBALLAGE / *PACKAGING*

Conditionnement :
Packaging :

Pochette plastique en polyéthylène (PE).
Polyethylene plastic pouch (PE).