

L A N D A N G E R	FICHE TECHNIQUE PRODUIT PRODUCT TECHNICAL SHEET	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-------------------	--	--

Editée le / Publishing date

2022-04-28

Page

1

REFERENCE	DESIGNATIONS - DESCRIPTIONS	CODE EAN
M03520	GABARIT DE CINTRAGE POUR PETITES PLAQUES 7 TROUS	3661219218784
	TEMPLATE BENDING FOR SMALL PLATES WITH 7 HOLES	

Catalogues / catalogs

Landanger 2015

Chapitre-page


14-020

DESCRIPTIF DU DISPOSITIF/ PRODUCT DESCRIPTION

L A N D A N G E R

Instrumentation de cambrage - Instruments for bending - Instrumentación para doblar

Fer à cintrer les mini plaques
Plate bender for mini plates
Llave para doblar mini placas

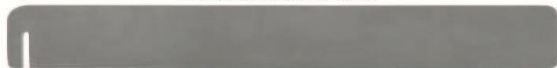

1,8 mm

M01250
12 cm - 4 3/4"

Fer à cintrer les plaques de reconstruction
Plate bender for reconstruction plates
Llave para doblar placas de reconstrucción



4,0 mm / 2,4 mm

M01251
14,5 cm - 5 3/4"

Fer à cintrer les plaques standards
Plate bender for standard plates
Llave para doblar placas rectas



5 mm

M01252
14,5 cm - 5 3/4"

4 mm / 2,4 mm

Fer à cintrer les plaques standards
Plate bender for standard plates
Llave para doblar placas rectas




5 mm

M01253
25 cm - 5 1/2"

7,3 mm

Gabarit de cintrage - Template bending - Modelo para doblar

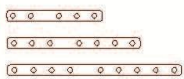

M01050
12 x 84 mm

Aluminium
Aluminum
Aluminio



M01050
12 x 84 mm

Petite plaque - Small plate - placa pequeña

M03520

M03530



XIV - Traumatologie
XIV - Traumatology
XIV Traumatología

Copyright © LANDANGER - 2009 - 2015

LANDANGER	FICHE TECHNIQUE PRODUIT <i>PRODUCT TECHNICAL SHEET</i>	FAQ 026 INFO Indice : 2 Date : 01/06/2018
-----------	---	---

Éditée le / *Publishing date*

2022-04-28

Page

2

REFERENCE	DESIGNATIONS - <i>DESCRIPTIONS</i>	CODE EAN
M03520	GABARIT DE CINTRAGE POUR PETITES PLAQUES 7 TROUS	3661219218784
	TEMPLATE BENDING FOR SMALL PLATES WITH 7 HOLES	

INFORMATIONS REGLEMENTAIRES / *REGULATORY INFORMATION*

Classe CE :
CE Class :

I

Code GMDN :
GMDN code :

44795

PERFORMANCES REVENDIQUEES / *INTENDED PERFORMANCE*

Dispositif médical de bloc opératoire :
Operating room medical device :

REUTILISABLE/REUSABLE

CARACTERISTIQUES SPECIFIQUE / *SPECIFIC FEATURES*

Type matière :
Material type :

ALUMINIUM

Nuance matière :
Material grade :

ALUMINIUM

Dimensions : Voir désignations produits
Outer dimensions : See products description
Chapitre-pages

Landanger 2015

14-020

Nettoyage/Pré-désinfection : Standard pour instrumentation de bloc opératoire réutilisable.
Cleaning / Pre-disinfection : Standard conditions for reusable surgical instruments

Stérilisation :

La stérilisation doit être effectuée selon un procédé de stérilisation à la vapeur conformément à la norme ISO 17665-1 (procédé sous vide fractionné) avec un stérilisateur conformément à la norme EN 285.

Méthode de stérilisation à la vapeur saturée	Température	Durée	Temps de séchage	Pression absolue
Pré-vide	132 °C (270 °F)	4 minutes	30 minutes	2,9 bars, 2,8 atm, 41,6 psi
Pré-vide	134 °C (273 °F)	3 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	5 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi
Pré-vide	134 °C (273 °F)	18 minutes	35 minutes	3,0 bars, 3,0 atm, 44,1 psi

Sterilization process :

The Sterilization shall be performed according to ISO 17665-1 validated steam sterilisation procedure (fractional vacuum procedure) in an autoclave in compliance with EN 285.

Saturated steam sterilization method	Temperature	Time	Drying time	Absolute pressure
Pre-vacuum	132 °C (270 °F)	4 minutes	30 minutes	2.9 bar, 2.8 atm, 41.6 psi
Pre-vacuum	134 °C (273 °F)	3 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	5 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi
Pre-vacuum	134 °C (273 °F)	18 minutes	35 minutes	3.0 bar, 3.0 atm, 44.1 psi

EMBALLAGE / *PACKAGING*

Conditionnement : Pochette plastique en polyéthylène (PE).
Packaging : Polyethylene plastic pouch (PE).