

Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

**GET YOURS HERE!**


Saisir le nom du produit, le numéro de lot, etc.

**777374** ▶ **Sigma-Aldrich®**

DBCO-Cy5

★★★★★ (0)

for Copper-free Click Chemistry

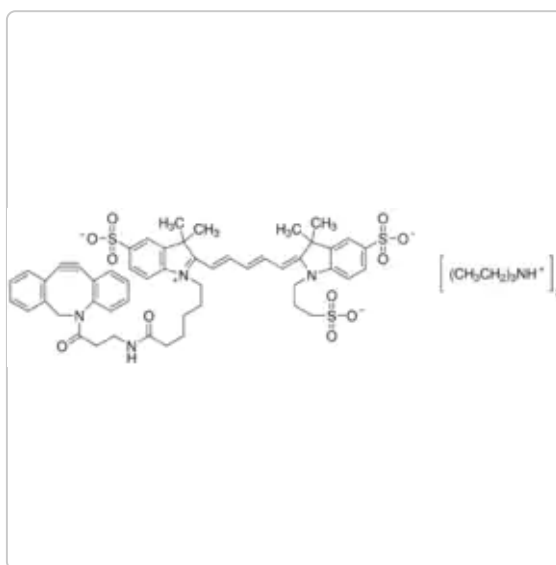[Toutes les photos](#) (3)**Synonyme(s):**

ADIBO-Cy5, Azadibenzocyclooctyne-Cy5, Dibenzocyclooctyne-Cy5

Empirical Formula (Hill Notation): $C_{64}H_{86}N_6O_{11}S_3$ **Poids moléculaire:**

1211.60

ID de substance PubChem:**329767903****NACRES:**

NA.22

HAUT

Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

GET YOURS HERE!



Reference

777374-1MG



1 MG

153,00 €



Disponible pour expédition le 10 octobre 2022 - **De**

—

1

+

Ajouter au panier

[Demander une commande en gros](#)

Documents


[↓ FDS](#)

[🔍 COO/COA](#)

[📄 Fiche des caractéristiques](#)

[Plus de documents >>](#)

PRODUITS RECOMMANDÉS

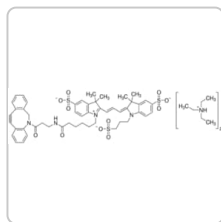

Sigma-Aldrich

777366

DBCO-Cy3

for Copper-free Click Chemistry

HAUT

[Consulter le prix et la disponibilité](#)

Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

GET YOURS HERE!

**PROPRIETES****Niveau de qualité**

100

Forme

solid

reaction suitability

reaction type: click chemistry

Temp. de stockage

-20°C

SMILES string
CC[NH+](CC)CC.O=C(CCNC(CCCCC[N+]=C(/C=C/C=C/C=C2C(C)(C)C3=CC(S([O-])
(=O)=O)=CC=C3N\2CCCS([O-])(=O)=O)C(C)(C)C4=C1C=CC(S([O-])
(=O)=O)=C4)=O)N5CC6=C(C=CC=C6)C#CC7=C5C=CC=C7.CC[NH+](CC)CC
InChI
1S/C52H56N4O11S3.2C6H15N/c1-51(2)42-34-40(69(62,63)64)25-27-45(42)54(31-14-6-9-22-
49(57)53-30-29-50(58)56-36-39-18-11-10-16-37(39)23-24-38-17-12-13-19-44(38)56)47(51)20-7-5-
8-21-48-52(3,4)43-35-41(70(65,66)67)26-28-46(43)55(48)32-15-33-68(59,60)61;2*1-4-7(5-2)6-
3/h5,7-8,10-13,16-21,25-28,34-35H,6,9,14-15,22,29-33,36H2,1-4H3,
(H3-,53,57,59,60,61,62,63,64,65,66,67);2*4-6H2,1-3H3
InChI key

PNVUJAWSQCGWCT-UHFFFAOYSA-N

Catégories apparentées[Click Chemistry Reagents](#)**DESCRIPTION**

HAUT
Application

This azadibenzocyclooctyne-cyanine dye derivative is a versatile labeling reagent for detection of azide containing molecules or compounds. Cyclooctynes are useful in strain-promoted copper-free azide-

Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

GET YOURS HERE!



Solubility : DMSO, DMF

Abs/Em - 646/661 nm

Extinction Coefficient:- 251,000 cm⁻¹M⁻¹

This azadibenzocyclooctyne-cyanine dye derivative is a versatile labeling reagent for detection of azide containing molecules or compounds. Cyclooctynes are useful in strain-promoted copper-free azide-alkyne cycloaddition reactions. This dibenzocyclooctyne will react with azide-functionalized compounds or biomolecules without the need for a Cu(I) catalyst to result in a stable triazole linkage.

Solubility : DMSO, DMF

Abs/Em - 646/661 nm

Extinction Coefficient:- 251,000 cm⁻¹M⁻¹

Conditionnement

1, 5 mg in amber glass bottle

INFORMATIONS SUR LA SÉCURITÉ

Pictograms



GHS07

Mention d'avertissement

Warning

Mentions de danger

H315 - H319 - H335

Conseils de prudence

P261 - P264 - P271 - P280 - P302 + P352 - P305 + P351 + P338

Classification des risques

Eye Irrit. 2 - Skin Irrit. 2 - STOT SE 3

Organes cibles

Respiratory system



Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

GET YOURS HERE!



WGK

WGK 3

Flash Point(F)

Not applicable

Point d'éclair C

Not applicable

DOCUMENTATION

Certificat d'analyse

Saisir un numéro de lot pour rechercher un certificat d'analyse (COA).

Numéro de lot

e.g. 023J5431

Comment saisir un numéro de lot (COA)

Rechercher

Certificat d'origine

Saisir un numéro de lot pour rechercher un certificat d'origine (COO).

Numéro de lot

e.g. 023J5431

Comment saisir un numéro de lot (COO)

Rechercher

HAUT

Plus De Documents

Structure Search


Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

GET YOURS HERE!



LES CLIENTS ONT ÉGALEMENT CONSULTÉ

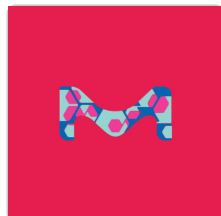

Sigma-Aldrich

GEPA15004

Cy[®]5 Bis NHS Ester

Cytiva PA15004, pack of 10 mg

[Consulter le prix et la disponibilité](#)

< 1 of 8 >

ARTICLES REVUS PAR DES PAIRS

Bio-orthogonal click-targeting nanocomposites for chemo-photothermal synergistic therapy in breast cancer.

Jianan Qiao et al.

Theranostics, 10(12), 5305-5321 (2020-05-07)

Chemo-photothermal synergistic treatment has a high potential to complement traditional cancer therapy and amplify its outcome. Precision in the delivery of these therapeutic agents to tumor cells has been indicated as being key to maximizing their therapeutic effects. Method: We

Targeting infected host cells in vivo via responsive azido-sugar mediated metabolic cell labeling followed by click reaction.

Yang Bo et al.

Biomaterials, 238, 119843-119843 (2020-02-18)

The early in vivo diagnosis of infectious disease foci is largely hindered by invasion and concealment of pathogens in host cells, making it difficult for conventional probes to detect and analyze intracellular pathogens. Taking advantage of the excessively produced reactive

Chemoenzymatic glycan labelling as a platform for site-specific IgM-antibody drug conjugate

Edward S X Moh et al.

Analytical biochemistry, 584, 113385-113385 (2019-08-06)

Immunoglobulin M (IgM) type antibodies play a significant role in complement activation, cellular debris clearance and cell quality control, and have the potential to be used as a therapeutic or targeting/delivery antibody. However, this potential has not been explored thoroughly

HAUT

Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

Molecular Cell, 80(5), 796-809 (2020-11-07)

The linkage, length, and architecture of ubiquitin (Ub) chains are all important variables in providing tight control over many biological paradigms. There are clear roles for branched architectures in regulating proteasome-mediated degradation, but the proteins that selectively recognize and process

GET YOURS HERE!



Replication Fork Activation Is Enabled by a Single-Stranded DNA Gate in CMG Helicase.

Michael R Wasserman et al.

Cell, 178(3), 600-611 (2019-07-28)

The eukaryotic replicative helicase CMG is a closed ring around double-stranded (ds)DNA at origins yet must transition to single-stranded (ss)DNA for helicase action. CMG must also handle repair intermediates, such as reversed forks that lack ssDNA. Here, using correlative single-molecule

[Afficher tous les articles scientifiques apparentés](#)


PROTOCOLES ET ARTICLES SCIENTIFIQUES

Articles

Copper-Free Click Chemistry

Copper-free click chemistry is an alternative approach to click chemistry that proceeds at a lower activation barrier and is free of cytotoxic transition metal catalysts.

PRODUITS CONSULTÉS RÉCEMMENT

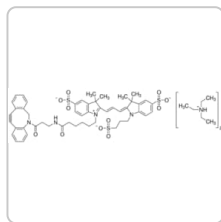

Sigma-Aldrich

777366

DBCO-Cy3

for Copper-free Click Chemistry

HAUT

[Consulter le prix et la disponibilité](#)

Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and

innovative products

GET YOURS HERE!



Soyez le premier à donner votre avis sur ce produit

Questions

Soyez le premier à poser une question

SERVICE TECHNIQUE

Notre équipe de scientifiques dispose d'une expérience dans tous les secteurs de la recherche, notamment en sciences de la vie, science des matériaux, synthèse chimique, chromatographie, analyse et dans de nombreux autres domaines..

Contactez notre Service technique

Le contenu de cette page vous a-t-il été utile ?*



Que pouvons-nous faire pour améliorer cette page spécifique de notre site Internet ?

Envoyer

HAUT

Scientific Newsletters

Subscribe to our eNewsletters & get tips on applications and innovative products

GET YOURS HERE!



La reproduction d'une quelconque partie du contenu de ce site est strictement interdite sans autorisation.

[Conditions d'utilisation du site](#) | [Politique de confidentialité](#) | [Conditions générales de vente](#) | [Consentement relatif au copyright](#) | [Paramètres des cookies](#)

