

How EZ Cytofunnels Work:

Designed for ease of use and efficiency, the single-use Thermo Scientific™ Shandon™ EZ Cytofunnel™ is an ideal solution for today's busy laboratories.

To use, insert the Thermo Scientific™ Shandon™ Cytoslide™ into the EZ Cytofunnel. Twelve EZ Cytofunnel assemblies fit into the Cytospin 4 head.




After cytocentrifugation, simply remove the EZ Cytofunnel from the Cytospin head. The assembly easily opens to remove the slide and dispose of the sample chamber.

A simple principle for excellent presentation

Thermo Scientific Cytofunnels and Slides

Thermo Scientific Shandon Cytofunnels

The original, trusted, easy-to-use funnel

	Description	Quantity	Catalog Number
	Thermo Scientific Shandon EZ Megafunnel	25 / Box	A78710001
	<ul style="list-style-type: none"> Designed for convenient preparation of larger volume samples, up to 6 ml Minimizes need to pre-centrifuge Cell deposition area of 22 X 14.75 mm (325 mm²) for increased cell yield Shandon™ Megafunnel Cytoslides™ included Can prepare both air-dried and fixed preparations providing valuable diagnostic information 		
	Thermo Scientific Shandon EZ Single Cytofunnel	40 / Box	A78710003
	With white filter cards, caps included	200 / Pack	A78710020
	Thermo Scientific Shandon EZ Single Cytofunnel	40 / Box	A78710004
	With brown filter cards, caps included		
	Thermo Scientific Shandon EZ Double Cytofunnel	40 / Box	A78710005
	With white filter cards, caps included	200 / Pack	A78710021
	<ul style="list-style-type: none"> Disposable sample chambers High-absorbency white filter cards can be used with sample volumes up to 0.5 ml Slower-absorbency brown filter cards can be used with sample volumes up to 0.4 ml Cell deposition area of 6 mm diameter (28 mm²) 		
	<ul style="list-style-type: none"> Disposable sample chambers Allows two sample deposition areas on a single slide Double the number of presentations per run Funnel markings correspond to slide markings for simple identification Pre-attached filter card 		


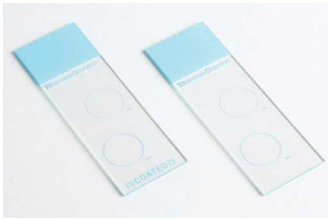
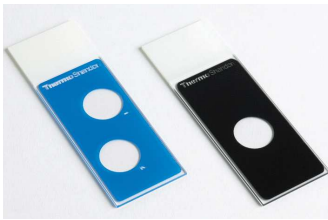
Thermo Scientific Shandon Cytoslides

Cytoslides

- Specially designed slides for the Cytospin 4
- Specimen area is outlined for rapid location during screening
- All slides are 75 x 25 mm, made from high-quality optical glass

Coated Slides

- Eliminates the need for albumin coating procedures
- Saves labor and resource time

	Description	Quantity	Catalog Number
	Thermo Scientific Shandon Single Cytoslides	Single circle, coated, circle on back	100 / Box 5991056
		Single circle, uncoated, circle on back	100 / Box 5991051
		Single circle, uncoated, circle on front	100 / Box 5991059
	<ul style="list-style-type: none"> Specifically designed for use with Shandon EZ Single Cytofunnels White 12 mm diameter circle on the slide defines where cells will be displayed to allow for efficient and rapid location 		
	Thermo Scientific Shandon Double Cytoslides	Double circle, coated	100 / Box 5991055
		Double circle, uncoated	100 / Box 5991054
	<ul style="list-style-type: none"> Specifically designed for use with Shandon EZ Double Cytofunnels Each specimen is displayed onto the slide within a circle marked "1" or "2" which corresponds to markings on the EZ Double Cytofunnel chamber 		
	Thermo Scientific Shandon Mask Cytoslides	Single circle, black mask	100 / Box 5991057
		Double circle, blue mask	100 / Box B5991050
	<ul style="list-style-type: none"> Specifically designed for immunofluorescence Black coating masks out UV light so that only the specimen within the circle can be viewed through a microscope 		