
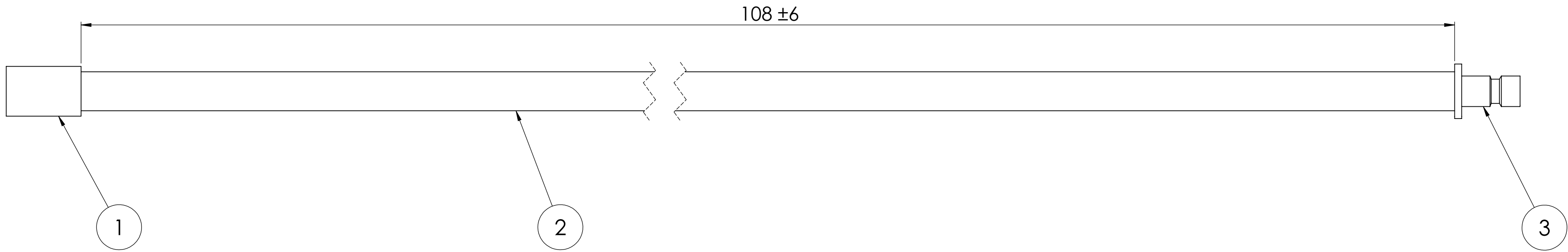


1REVISIONS		
REV	ECO NUMBER	NAME
A	PR056827	M. GROBMAN

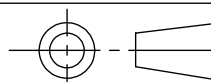

ASSEMBLY INSTRUCTIONS:

1. INSERT BARBED END OF FITTINGS INTO THE HOSE. HOSE ENDS FIT FLUSH TO END OF FITTINGS.
2. THE AIR LEAKAGE SHALL NOT EXCEED A PRESSURE DROP OF 4 mmHg AT 300 mmHg FOR 1 MINUTE.
3. NEATLY COIL PRODUCT FOR PACKAGING.

THIS DOCUMENT CONTAINS PROPRIETARY TRADE SECRET AND CONFIDENTIAL INFORMATION WHICH IS THE PROPERTY OF SPACELABS HEALTHCARE. THIS DOCUMENT AND THE INFORMATION IT CONTAINS ARE NOT TO BE COPIED, DISTRIBUTED, OR DISCLOSED TO OTHERS, OR USED IN ANY MANNER OUTSIDE OF SPACELABS HEALTHCARE WITHOUT PRIOR WRITTEN APPROVAL OF SPACELABS HEALTHCARE.			Spacelabs Healthcare		
UNLESS OTHERWISE SPECIFIED DIMENSIONS ARE IN INCHES. TOLERANCES AS SHOWN BELOW			DESCRIPTION		
LINEAR ANGULAR X.X ± 0.25 ± 1 ° X.XX ± 0.25 X.XXX ± 0.125 X.XXXX ± 0.0625			INTERPRET DIMENSIONS AND TOLERANCES PER ISO OR ASME Y14.5M-1994 DIMENSIONS SHOWN ARE FINISHED SIZES AND SHALL ALLOW FOR PLATING AND COATING, BUT NOT PAINTING		
 THIRD ANGLE PROJECTION			GENERAL THREAD FIT 6H 6g		
SIZE		PART NUMBER		REV	
A1		653-1660-00		A	
VERSION:				SHEET 1 OF 1	